

Start-up Package for Single Station Personal Pipettor

Xpp-single

All you need to start automated pipetting of samples into 96 or 384-well microplates.

With a **15% discount** on all packages, prices start at

£12,610 / €18,280



delivering solutions

- Each Xpp-Single Start-up Package includes:
 - A Xpp-Single high performance automated multichannel pipettor.
Three models to choose from: **Xpp-Single96L** (low volume 96 channel),
Xpp-Single96H (high volume 96 channel),
Xpp-Single384 (384 channel version).
 - Corresponding EZ-Load pipetting head.
 - 2 cases of EZ-Load pipette tips.
- The Xpp-Single automated multichannel pipettor uses a microprocessor-controlled single plate station which provides accurate and precise pipetting of microliter volumes of samples into 96 or 384-well microplates.
- Small footprint conserves bench space and fits easily into laminar flow and class II cabinets.
- Compatible with EZ-Load pipettor heads to simplify the use of disposable tips for multiple channel pipetting. Tip loading and unloading becomes a simple manual process.
- Suitable for plate replication, plate reformatting, reagent addition, compound addition and serial dilution (by row or column).
- No PC Required. Simple to use on-board programming with multi-line display.
- **EZ-Load Pipetting Heads**
Specifically-designed tip racks enable a rack of tips to be transferred directly from the package into the EZ-Load head. Simply drop the rack of tips into the head and place the head in the pipettor.
- A reliable seal is provided for optimal pipetting performance. The top of the pipette tips are sealed against the face seal by the pipettor's automated loading system. Compatible with a wide range of other 96 and 384 channel pipettors, including all other fluidX Xpp and Xps systems.
- Capable of nanolitre dispensing when used with FluidX PocketTips!

Start-up Packages for Single Station Personal Pipettor

XPP SINGLE SPECIFICATIONS

Dispensing Precision	<2% CV at 1 ul (Xpp-Single96L), <2% CV at 5 ul (Xpp-Single96H) <2.8%CV at 1 ul (Xpp-Single384)
Dispensing Accuracy	+/- 1% at 1 ul (Xpp-Single96L), +/- 1% at 10 ul (Xpp-Single96H) +/- 1.25% at 1 ul (Xpp-Single384)
Resolution	1 ul (Xpp-Single96H) 0.1 ul (Xpp-Single96L & Xpp-Single384)
Dimensions	12"L x 11"W x 23"H 304L x 280W x 585H mm
Weight Approx	40 lbs (18.2 kg)
Power	100 to 240 Vac, 50/60 Hz
Operating Temperature	10 - 38°C

PACKAGE 1 - Low Volume 96-Channel System

£12,610 / €18,280 (15% discount on list price)

Cat. No.	Description
30-1900	Xpp-Single96L, 96-channel low-volume version with 125ul pipetting channels for optimal pipetting between 0.5ul and 125ul. Motorised tip locking mechanism.
30-9002	"EZ Load" pipetting head, 96-channel, for fluidx ESP tips.
36-1001	96 EZ-Load pipette tips, 125ul, non-sterile, 96 tips per cartridge, 10 cartridges per case.

PACKAGE 2 - High Volume 96-Channel System

£12,800 / €18,550 (15% discount on list price)

Cat. No.	Description
30-1901	Xpp-Single96H, high-volume 96-channel version with 550ul pipetting channels for optimal pipetting between 5ul and 550ul. Motorised tip locking mechanism.
30-9001	"EZ Load" Pipetting Head, 96-Channel, for fluidx 550ul tips.
35-1001	96 EZ-Load Pipette Tips, 500ul, non sterile, 96 tips per cartridge, 20 cartridges per case.

PACKAGE 3 - 384-Channel System

£20,970 / €30,410 (15% discount on list price)

Cat. No.	Description
30-1902	Xpp-Single384, 384-channel version with for full 384 plate dispensing and replication and additional 96-channel head compatibility. Motorised tip locking mechanism.
30-9004	"EZ Load" pipetting head, 384-channel, for 384 fluidx ESP Tips.
36-2001	384 EZ-Load pipette tips, 125ul, Non-sterile. 10 racks per case.


For further information, or to place an order simply contact FluidX on

Tel: +44 (0) 1625 861614

Email: order@fluidx.co.uk

Standard terms and conditions apply.

FluidX reserve the right to update and change product information without notice.

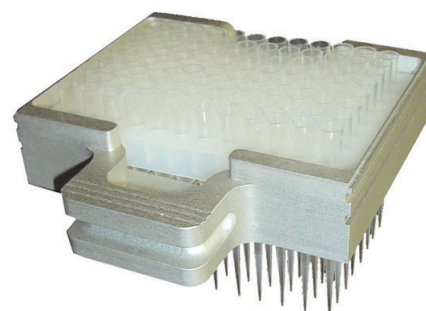

Tel:+44 (0) 1625 861614 Email:info@fluidx.co.uk www.fluidx.co.uk